

FlexiGuide Needles

Sterile / Disposable

Catalog # 0911 / 0912 / 0919



FlexiGuide Needles used in HDR Brachytherapy:

Prostate • Interstitial / Intracavitary Cervix • Breast • Others

Provides Smooth and Accurate Needle Placement within the Gland,
Conforms to the Shape of the Prostate Gland,
Minimal Trauma and given its Flexibility provides Minimal Discomfort for Patient,
CT Compatibility provides excellent Delineation for Treatment Planning Contouring with Minimal CT Artifacts.

GAUGE	LENGTH	CAT. #
17	25 cm	0919
15.5	20 cm	0911
15.5	25 cm	0912

Features:

- FlexiGuide Needles are Sterile & Disposable
- Ripple-Hubs for better Grip
- Color-coded for Size and Length Identification
- Rigid Obturators for Needle Stabilization
- Pencil-point Sharpened
- Configured for Varian's® Click-fit Connections

Packaged:

25 Needles / Box; 3 Pouches - 5 Needles ea; 10 Pouches - 1 Needle each.

Compatibility:

GammaMed® and VariSource® Remote Afterloading Systems

An ISO 13485:2003 Certified Company



Mick Radio-Nuclear Instruments, Inc.

"The Brachytherapy Company"


521 Homestead Avenue, Mt. Vernon, NY, 10550, U.S.A.

Tel: 914-667-3999 • Fax: 914-665-8834 • Toll-Free: 877-597-MRNI

E-mail: sales@micknuclear.com • Website: www.micknuclear.com


HDR Templates used in Conjunction with FlexiGuide Needles

HDR Contour Template - Small Matrix *"PROSTATE" 15.5 Gauge FlexiGuide Needles*

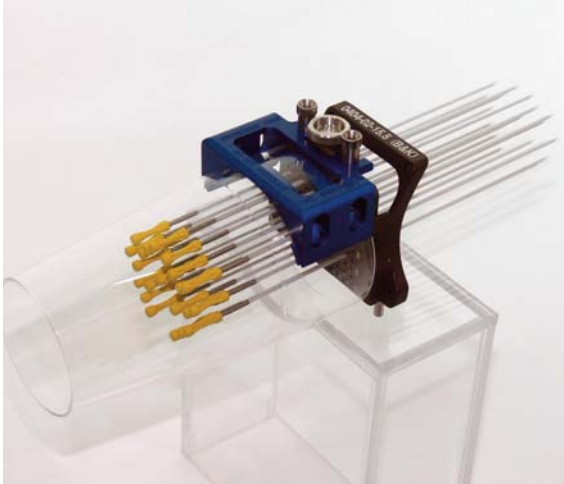

Small matrix HDR Contour Template shown with 15.5 Gauge FlexiGuide Needles (Yellow), Rigid "Ripple-Hub" Stylets and attached Protective Acrylic Cone.
B&K Compatible.

HDR Contour Template - Large Matrix *"PROSTATE" 17 Gauge FlexiGuide Needles*

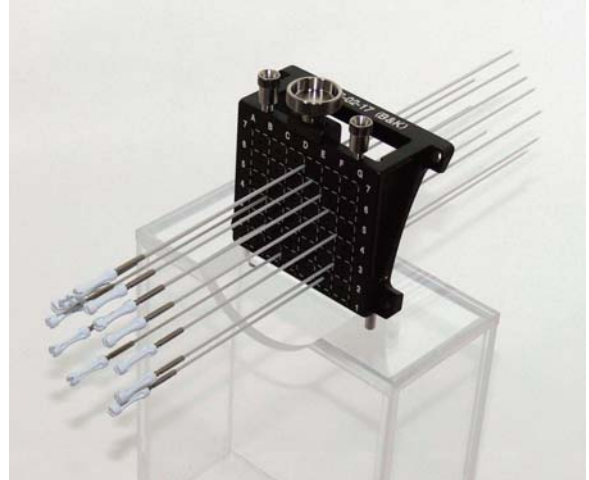


Large matrix HDR Contour Template shown with 17 Gauge FlexiGuide Needles (White), Rigid "Ripple-Hub" Stylets and without attached Protective Acrylic Cone.
B&K Compatible.

HDR MAC Template - Small Matrix *"Vaginal / Cervical" 15.5 Gauge FlexiGuide Needles*



Small matrix HDR MAC Template shown with Vaginal Obturator, Cervical Tandem and laterally disbursed 20 cm FlexiGuide Needles.
The Tandem is replaceable with a 25 cm long FlexiGuide Needle (Green).

HDR MAC Template - Large Matrix *"Vaginal / Cervical" 15.5 Gauge FlexiGuide Needles*

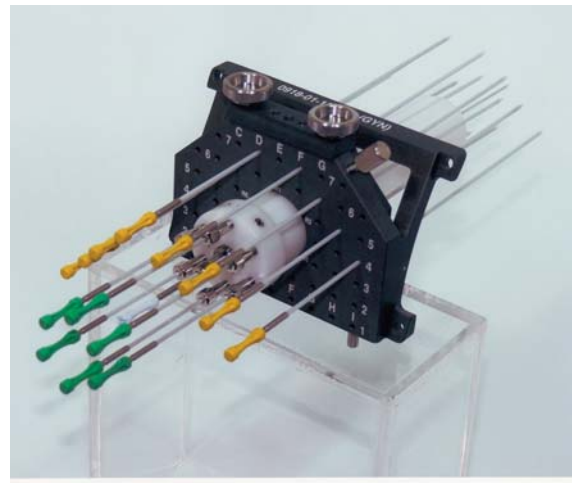

Large matrix HDR MAC Template shown with Vaginal Obturator, Cervical Tandem and "extended" laterally disbursed 20 cm FlexiGuide Needles.
The Tandem is replaceable with a 25 cm long FlexiGuide Needle (Green) as shown.

For more information, please visit our website;
www.micknuclear.com
or call 914-667-3999