


CT/MR Contour Prostate Template Sets

Sets 1204-02-XX, 1204-04-XX, 1206-02-XX, 1206-04-XX

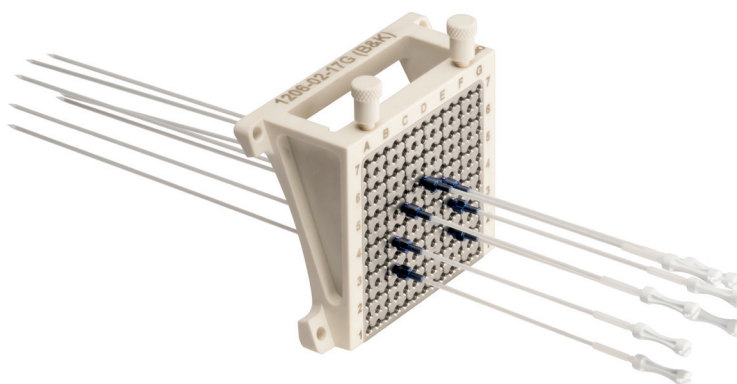

Facilitates interstitial needles for the treatment of the prostate.

Prostate

Simplicity

The CT/MR Contour Prostate Template facilitates interstitial needles for the treatment of the prostate. All templates are laser engraved with letters to denote the different vertical columns of needle channels and with numbers to denote the different horizontal rows of needle channels. The numbering/lettering layout corresponds to the type of stepper being used.

Flexibility

The CT/MR Contour Prostate Template is available to accommodate needles of different gauges which are securely held in place by color coded needle collets associated with the appropriate needle size. Small and large template sizes are provided. The CT/MR Contour Prostate Template is curved to fit the sloping anatomy of the perineum. There are four suture holes provided for optional fixation of the template to the perineum.


Reliability

The CT/MR Contour Prostate Template is constructed of PPSU (Polyphenylsulfone) and the template inserts and needle collets are constructed of titanium, which is ideal for continuous use and steam sterilization. Furthermore, protective cones are available for the small and the large templates to shield the implanted needles from accidental displacement or impact.

Features

- Contoured to fit the sloping perineum
- Reusable and autoclavable
- MR conditionality allows for 1.5T or 3T

- Individual locking collet design enhances discrete needle targeting
- Available for variety of needle sizes: 15.5G, 17G, 18G, 1.65 mm, 2 mm
- Provided with four suture holes allowing the template to be directly sutured to the perineum
- Resistance-free bore holes for smooth needle entry
- Durable construction
- Optional Protective Cone

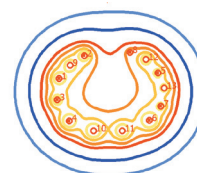

Isodose



Template Needle Collets



Eckert & Ziegler

Mick Radio-Nuclear
Instruments, Inc.

An Eckert & Ziegler BEBIG Company

CT/MR Contour Prostate Template Sets

Set 1204-02-XX (Small Matrix), for CIVCO Steppers, BK labeling scheme¹

Qty.	Item	Part number
1	CT/MR Contour Prostate Template PPSU Body, small matrix	
1	CT/MR Contour Prostate Template Titanium Insert, small matrix	
1	Titanium Collet Wrench	1205-02
1	Mounting Pegs, (2 pcs)	1205-03
1	Mounting Bracket Lock Knobs (2 pcs) for Protective Cone, reusable	1205-12
Choose one CT/MR Contour Prostate Template, small matrix, BK based on needle size (XX)		
1	CT/MR Contour Prostate Template, small matrix, BK, 15.5G	1204-02-15.5
1	CT/MR Contour Prostate Template, small matrix, BK, 17G	1204-02-17
1	CT/MR Contour Prostate Template, small matrix, BK, 17G, Eckert & Ziegler BEBIG	1204-02-EZ-17
1	CT/MR Contour Prostate Template, small matrix, BK, 18G	1204-02-18
1	CT/MR Contour Prostate Template, small matrix, BK, 1.65 mm	1204-02-1.65
1	CT/MR Contour Prostate Template, small matrix, BK, 1.65 mm, Eckert & Ziegler BEBIG	1204-02-EZ-1.65
1	CT/MR Contour Prostate Template, small matrix, BK, 2.0 mm	1204-02-2.0
Choose one Template Needle Collets Set (36 pcs) based on needle size (XX)		
6	Template Needle Collets (Titanium), gray, 15.5G, (6 pcs)	1105-01-15.5
6	Template Needle Collets (Titanium), blue, 17G, (6 pcs)	1105-01-17
6	Template Needle Collets (Titanium), seafoam, 17G, Eckert & Ziegler BEBIG, (6 pcs)	1105-01-EZ-17
6	Template Needle Collets (Titanium), gold, 18G, (6 pcs)	1105-01-18
6	Template Needle Collets (Titanium), magenta, 1.65 mm, (6 pcs)	1105-01-1.65
6	Template Needle Collets (Titanium), bronze, 1.65 mm, Eckert & Ziegler BEBIG, (6 pcs)	1105-01-EZ-1.65
6	Template Needle Collets (Titanium), green, 2.0 mm, (6 pcs)	1105-01-2.0

Set 1204-04-XX (Small Matrix), for CIVCO Steppers, Siemens labeling scheme¹

Qty.	Item	Part number
1	CT/MR Contour Prostate Template PPSU Body, small matrix	
1	CT/MR Contour Prostate Template Titanium Insert, small matrix	
1	Titanium Collet Wrench	1205-02
1	Mounting Pegs, (2 pcs)	1205-03
1	Mounting Bracket Lock Knobs (2 pcs) for Protective Cone, reusable	1205-12
Choose one CT/MR Contour Prostate Template, small matrix, Siemens based on needle size (XX)		
1	CT/MR Contour Prostate Template, small matrix, Siemens, 15.5G	1204-04-15.5
1	CT/MR Contour Prostate Template, small matrix, Siemens, 17G	1204-04-17
1	CT/MR Contour Prostate Template, small matrix, Siemens, 17G, Eckert & Ziegler BEBIG	1204-04-EZ-17
1	CT/MR Contour Prostate Template, small matrix, Siemens, 18G	1204-04-18
1	CT/MR Contour Prostate Template, small matrix, Siemens, 1.65 mm	1204-04-1.65
1	CT/MR Contour Prostate Template, small matrix, Siemens, 1.65 mm, Eckert & Ziegler BEBIG	1204-04-EZ-1.65
1	CT/MR Contour Prostate Template, small matrix, Siemens, 2.0 mm	1204-04-2.0
Choose one Template Needle Collets Set (36 pcs) based on needle size (XX)		
6	Template Needle Collets (Titanium), gray, 15.5G, (6 pcs)	1105-01-15.5
6	Template Needle Collets (Titanium), blue, 17G, (6 pcs)	1105-01-17
6	Template Needle Collets (Titanium), seafoam, 17G, Eckert & Ziegler BEBIG, (6 pcs)	1105-01-EZ-17
6	Template Needle Collets (Titanium), gold, 18G, (6 pcs)	1105-01-18
6	Template Needle Collets (Titanium), magenta, 1.65 mm, (6 pcs)	1105-01-1.65
6	Template Needle Collets (Titanium), bronze, 1.65 mm, Eckert & Ziegler BEBIG, (6 pcs)	1105-01-EZ-1.65
6	Template Needle Collets (Titanium), green, 2.0 mm, (6 pcs)	1105-01-2.0

CT/MR Contour Prostate Template Sets

Set 1206-02-XX (Large Matrix), for CIVCO Steppers, BK labeling scheme¹

Qty.	Item	Part number
1	CT/MR Contour Prostate Template PPSU Body, large matrix	
1	CT/MR Contour Prostate Template Titanium Insert, large matrix	
1	Titanium Collet Wrench	1205-02
1	Mounting Pegs, (2 pcs)	1205-03
1	Mounting Bracket Lock Knobs (2 pcs) for Protective Cone, reusable	1205-12
Choose one CT/MR Contour Prostate Template, large matrix, BK based on needle size (XX)		
1	CT/MR Contour Prostate Template, large matrix, BK, 15.5G	1206-02-15.5
1	CT/MR Contour Prostate Template, large matrix, BK, 17G	1206-02-17
1	CT/MR Contour Prostate Template, large matrix, BK, 17G, Eckert & Ziegler BEBIG	1206-02-EZ-17
1	CT/MR Contour Prostate Template, large matrix, BK, 18G	1206-02-18
1	CT/MR Contour Prostate Template, large matrix, BK, 1.65 mm	1206-02-1.65
1	CT/MR Contour Prostate Template, large matrix, BK, 1.65 mm, Eckert & Ziegler BEBIG	1206-02-EZ-1.65
1	CT/MR Contour Prostate Template, large matrix, BK, 2.0 mm	1206-02-2.0
Choose one Template Needle Collets Set (36 pcs) based on needle size (XX)		
6	Template Needle Collets (Titanium), gray, 15.5G, (6 pcs)	1105-01-15.5
6	Template Needle Collets (Titanium), blue, 17G, (6 pcs)	1105-01-17
6	Template Needle Collets (Titanium), seafoam, 17G, Eckert & Ziegler BEBIG, (6 pcs)	1105-01-EZ-17
6	Template Needle Collets (Titanium), gold, 18G, (6 pcs)	1105-01-18
6	Template Needle Collets (Titanium), magenta, 1.65 mm, (6 pcs)	1105-01-1.65
6	Template Needle Collets (Titanium), bronze, 1.65 mm, Eckert & Ziegler BEBIG, (6 pcs)	1105-01-EZ-1.65
6	Template Needle Collets (Titanium), green, 2.0 mm, (6 pcs)	1105-01-2.0

Set 1206-04-XX (Large Matrix), for CIVCO Steppers, Siemens labeling scheme¹

Qty.	Item	Part number
1	CT/MR Contour Prostate Template PPSU Body, large matrix	
1	CT/MR Contour Prostate Template Titanium Insert, large matrix	
1	Titanium Collet Wrench	1205-02
1	Mounting Pegs, (2 pcs)	1205-03
1	Mounting Bracket Lock Knobs (2 pcs) for Protective Cone, reusable	1205-12
Choose one CT/MR Contour Prostate Template, large matrix, Siemens based on needle size (XX)		
1	CT/MR Contour Prostate Template, large matrix, Siemens, 15.5G	1206-04-15.5
1	CT/MR Contour Prostate Template, large matrix, Siemens, 17G	1206-04-17
1	CT/MR Contour Prostate Template, large matrix, Siemens, 17G, Eckert & Ziegler BEBIG	1206-04-EZ-17
1	CT/MR Contour Prostate Template, large matrix, Siemens, 18G	1206-04-18
1	CT/MR Contour Prostate Template, large matrix, Siemens, 1.65 mm	1206-04-1.65
1	CT/MR Contour Prostate Template, large matrix, Siemens, 1.65 mm, Eckert & Ziegler BEBIG	1206-04-EZ-1.65
1	CT/MR Contour Prostate Template, large matrix, Siemens, 2.0 mm	1206-04-2.0
Choose one Template Needle Collets Set (36 pcs) based on needle size (XX)		
6	Template Needle Collets (Titanium), gray, 15.5G, (6 pcs)	1105-01-15.5
6	Template Needle Collets (Titanium), blue, 17G, (6 pcs)	1105-01-17
6	Template Needle Collets (Titanium), seafoam, 17G, Eckert & Ziegler BEBIG, (6 pcs)	1105-01-EZ-17
6	Template Needle Collets (Titanium), gold, 18G, (6 pcs)	1105-01-18
6	Template Needle Collets (Titanium), magenta, 1.65 mm, (6 pcs)	1105-01-1.65
6	Template Needle Collets (Titanium), bronze, 1.65 mm, Eckert & Ziegler BEBIG, (6 pcs)	1105-01-EZ-1.65
6	Template Needle Collets (Titanium), green, 2.0 mm, (6 pcs)	1105-01-2.0

CT/MR Contour Prostate Template Sets

Required Items Sets 1204-02-EZ-17, 1204-04-EZ-17, 1206-02-EZ-17, 1206-04-EZ-17

Qty.	Item	Part number
1	Disposable Steel Needle Set	SET0052
20	Transfer Tube, D = 3 mm, L = 1000 mm	LAF1000

Required Items Sets 1204-02-EZ-1.65, 1204-04-EZ-1.65, 1206-02-EZ-1.65, 1206-04-EZ-1.65

Qty.	Item	Part number
1	Plastic Needle Set	SET0012
20	Easy Click Transfer Tube, D = 3 mm, L = 1000 mm	LAG1000



Eckert & Ziegler BEBIG SagiNova® afterloader

Optional Items

Qty.	Item	Part number
1	Protective Cone with Mounting Bracket (Small, Aluminium), reusable	9807-PC II
1	Protective Cone with Mounting Bracket (Small, PPSU), reusable	1204-PC
1	Protective Cone with Mounting Bracket (Large, Aluminium), reusable	9807-PC
1	Protective Cone with Mounting Bracket (Large, PPSU), reusable	1206-PC
1	Mounting Bracket for Protective Cone, (Small, Aluminium), reusable	1205-11-01
1	Mounting Bracket for Protective Cone, (Small, PPSU), reusable	1204-04
1	Mounting Bracket for Protective Cone, (Large, Aluminium), reusable	1205-11-02
1	Mounting Bracket for Protective Cone, (Large, PPSU), reusable	1206-04
1	Protective Cone, reusable	9807-04
1	Basket with lid for cleaning and sterilizing of small parts, 121 x 73 x 36 mm	LLZ07-01
1	Sterilization Box, 381 x 254 x 38 mm	96-1522A

- ¹ The CT/MR Contour Prostate Template has been tested for compatibility with the following accessories and external devices based on stated nominal dimensions. Users should validate that their particular needle is compatible with the template and collet prior to use:
- Plastic & metal needles diameters that can be used with the corresponding template include: 15.5G, 17G, 18G, 1.65 mm (16G) and 2.0 mm

Federal (USA) law restricts the sale of this device by or on order of a physician. For a complete listing of all indications, contraindications, precautions and warnings, please refer to the Instruction for Use which accompany each product.

Mick and SagiNova are registered trademarks of Eckert & Ziegler BEBIG SA and its subsidiaries.

The mentioned products are not available in all markets. Please contact your local Eckert & Ziegler BEBIG representative for more information.

Corporate Head Office:

**Eckert & Ziegler
BEBIG SA**
Rue Jules Bordet
Zone Industrielle C
7180 Seneffe
Belgium

Phone +32 64 520 811
Fax +32 64 520 801
info@bebig.com

Manufacturer:

**Mick Radio-Nuclear Instruments, Inc.
An Eckert & Ziegler BEBIG Company**
521 Homestead Avenue
Mount Vernon, NY 10550
USA

Phone +1914 667 3999
Fax +1914 665 8834
sales@micknuclear.com

Regional Sales, Marketing and Service:

Europe, Middle East, Africa,
Latin America, Asia Pacific

**Eckert & Ziegler
BEBIG SA**
Rue Jules Bordet
Zone Industrielle C
7180 Seneffe
Belgium

Phone +32 64 520 811
Fax +32 64 520 801
info@bebig.com

North America

**Mick Radio-Nuclear Instruments, Inc.
An Eckert & Ziegler BEBIG Company**
521 Homestead Avenue
Mount Vernon, NY 10550
USA

Phone +1914 667 3999
Fax +1914 665 8834
sales@micknuclear.com

www.bebig.com
www.micknuclear.com